

American Health Gate

DISPOSABLE STERILE SAFETY LANCETS AHG-IIA

Instruction Manual



TABLE OF CONTENTS

Product Description	3
Intended Purpose	3
Indications for Use	3
Intended User & Patient Population	4
Contraindications	4
Warnings and Precautions	4
Instructions for Use	5
Storage Conditions	6
Shelf Life	6
Index of Symbols	7

PRODUCT DESCRIPTION

AHG- IIA DISPOSABLE STERILE SAFETY LANCETS

The AHG Disposable Sterile Safety Lancet (Model: AHG-IIA) is a sterile, single-use medical device intended for capillary blood sampling. It features a pre-loaded needle that automatically retracts after activation to reduce the risk of needlestick injury and prevent reuse.

Specifications

- Gauge: 28G
- Penetration depth: 1.8 mm
- Needle diameter: 0.36 mm
- Sterilization method: Irradiation

INTENDED PURPOSE

The AHG-IIA is intended to puncture the skin (e.g., fingertip or other appropriate capillary sampling site) to obtain a small volume of capillary whole blood for in vitro diagnostic testing (e.g., blood glucose monitoring).

The device is intended for single-patient use only.

INDICATIONS FOR USE

Indicated for the collection of capillary whole blood samples for diagnostic testing, including monitoring of physiological parameters such as blood glucose levels.

INTENDED USER & PATIENT POPULATION

The device is intended for:

- Healthcare professionals in clinical or laboratory settings.
- Lay users performing self-testing at home.

CONTRAINDICATIONS


- Do not use in areas with infection, edema, or compromised skin integrity.
- Use with caution in patients with bleeding disorders or on anticoagulant therapy.
- Not intended for use on infants unless directed by a healthcare professional.

WARNINGS AND PRECAUTIONS

- For single use only. Do not reuse or resterilize.
- Do not use after the expiry date indicated on the packaging.
- Do not use if the protective cap is missing or damaged, or if sterile barrier packaging is opened, damaged, or compromised.
- Wash hands thoroughly with soap and warm water before use.
- Use a new lancet for each puncture.
- Used lancets are potentially biohazardous.
- Improper use may increase the risk of infection or injury.
- Keep out of reach of children.
- In the event of a serious incident, report it to the manufacturer and the competent authority in your EU Member State.


INSTRUCTIONS FOR USE

Note: Wash your hands with soap and warm water and dry thoroughly before and after use.

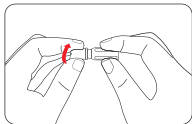


Step 1 – Preparation

Twist off the protective cap.
Do not pull straight off.

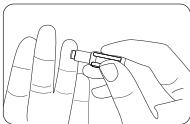


Step 2 – Position

Place the lancet opening firmly
against the puncture site.

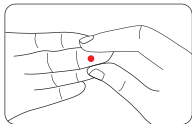

Step 3 – Activation

Press the device down firmly.
The needle will deploy and
retract automatically.

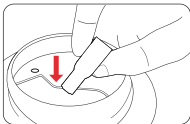

Step 4 – Disposal

Immediately dispose of the
used device in an approved
sharps container.

▲ Potential Biohazard

Used lancets may be contaminated with blood. Handle and dispose of in accordance with local regulations for sharps and biohazardous waste.

STORAGE CONDITIONS

- Store in a cool, dry place.
- Do not expose to temperatures above 30°C (86°F).
- Protect from moisture and direct sunlight.

SHELF LIFE

Shelf life: 5 years from date of manufacture, if stored under recommended conditions. Do not use after the expiration date.

DISPOSAL

Dispose of used lancets as medical sharps and biohazardous waste in accordance with local/national regulations.

INDEX OF SYMBOLS

The following symbols may appear on the packaging and in these instructions for use.

	Lot Number		Consult Instructions for Use
	Catalogue Number		Do Not Reuse
	Date of Manufacturer		Manufacturer
	Use-by Date		Sterilized Using Irradiation
	Keep Away from Sunlight		Medical Device
	Keep Dry		Do Not Resterilize
	Caution		Importer
	CE Marking		Unique Device Identification
	Authorized Representative in the EU		
	Dispose of in accordance with local regulations		



Head office:

American Health Gate LLC

450 Alaskan Way S., Suite 200, Seattle, WA 98104, USA

Website: www.americanhealthgate.com